

KRYPTON AND XENON ANALYSIS

Technique : Micro-GC

Column : Tamis


Carrier gas : Helium

Injector : BF

Column temperature : 120°C

Column pressure : 28 psi

μTCD sensitivity : Standard

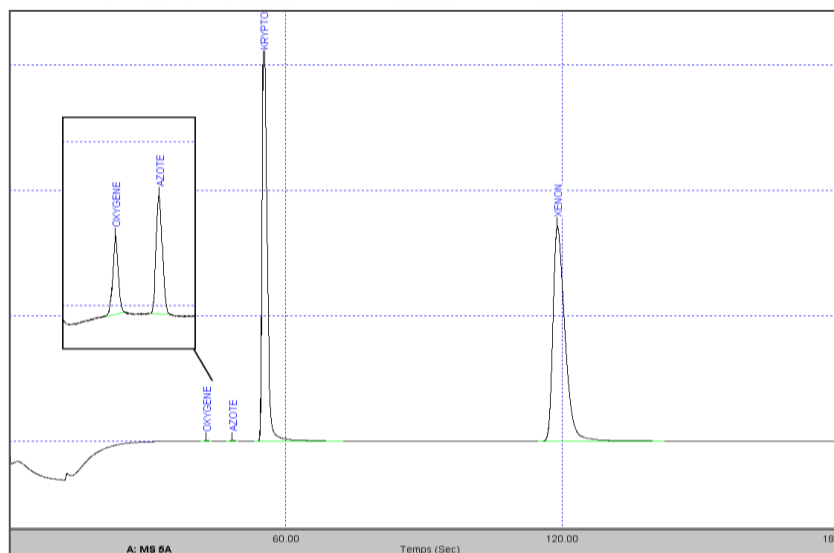

The Micro-GC is a powerful solution that provides fast, accurate and easy-to-use solution to analyze your gas sample.

A wide range of capillary columns and MEMS injectors are available to optimize channel combinations inside the Micro-GC to obtain the desired separation. The μ TCD is the high-performance universal detector used in such systems. A parallel analytical configuration of the channels or analytical modules will allow the sampling of a small amount of gas sample and simultaneous injection in all channels to develop up to 4 chromatograms and one cumulative report..

SRA Instruments is the ideal Partner with specific competences to develop with you turn key solutions based on optimized Micro-GC, accessories and dedicated software utility or report calculations for specific plant context.

Some examples of applications are:
gas coke analysis,
biogas, syngas,
refinery and fuel gas,
natural gas, fuel cells.



More information : www.sra-instruments.com